


Gen 5 / Optima Driver Installation

1. Make sure the amplifier and force plate are connected and the amplifier is powered on.
2. In the control panel search for the the '**Device Manager**' and click on it.
3. Locate the Gen 5 device. It will likely be listed under '**Other devices**' and have a yellow warning on its logo.
4. Right click the GEN 5 device and select '**Update Driver Software ...**'

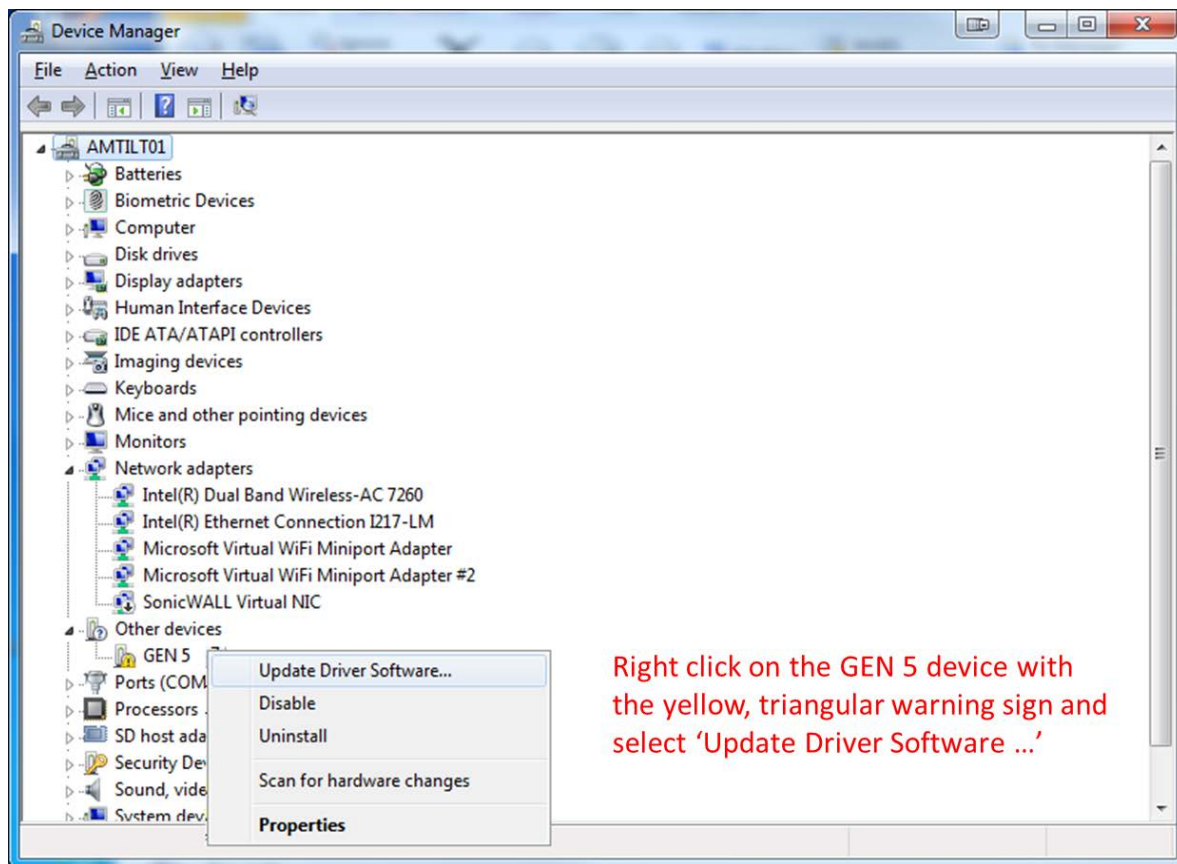

5. Select '**Browse my computer for driver software**'.
6. Click '**Browse**' and navigate to the folder containing the driver files. Make sure to select the appropriate 32-bit or 64-bit version for your operating system and then click '**OK**'.
 - a. To determine which operating system you have click '**Computer**' under the Windows '**Start**' menu, then right click a select '**Properties**'.
7. Click '**OK**' and then click '**Next**'
8. The drivers should install.

For more detailed instructions on the driver installation see page 17 of the Gen 5 user manual.